

H2 Series Product's
Professionally Manufactured



OPERATING INSTRUCTION

Hydrogen Eye Massager (User manual)

MR-03



Manniskans Resurser AB
Stockholm, Sweden
Tel: +468334908
mail@mr-ab.se
humans-resources.one
aranovich.se

This product is not a medical device and cannot
be advertised for any therapeutic effect.




Warranty Card

Name: _____ Mobile Contact: _____

Mailing Address: _____


Fault Description _____

④

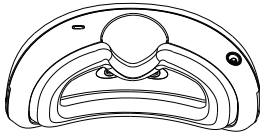

After each use, pour out the condensed water in the mist hole and dry it with a paper towel.

⑤

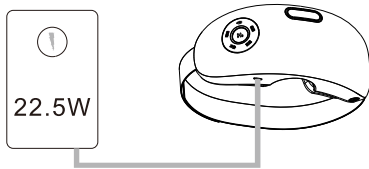

Use the "hot compress function", when the amount of fog is in sufficient, please fully charge it before use; if you need to use a power bank, it is recommended to use a 22.5W charging power supply.

Tips

- This product is not a medical device.
- Keep the liquid tank clean(it is recommended to clean it with pure water before and after each use).
- For specific operations, please follow the instructions in the manual.

Catalogue	1
Product Description	2-4
Safety Precautions	4-6
Parts list	6
Product Parameters	7
Structure and Panel Introduction	8-9
How to use	9-15
Maintenance and Storage	16-18
Troubleshooting	19-20
After-sales Service	21-22
Operation Process Card	23-24

respected user:

Sincerely thank you for choosing this product. The "Instruction Manual" will help you master the use of this product as soon as possible, allowing you to fully experience the enjoyment of a high-quality life. Before use, please read this manual carefully and pay special attention to safety matters to ensure correct operation and use. If you encounter any problems during use, please contact our company's after-sales service center directly. We will serve you wholeheartedly.

Product Description

Product Usage:


1. It is suitable for alleviating or reducing visual fatigue caused by overuse of eyes such as work, study, driving, staying up late, and environmental changes.
2. Moisturizing, hot compress and cold compress physical therapy for dry eyes, pain and other uncomfortable symptoms caused by dry eye.
3. Promote blood circulation around the eyes, reduce dark circles and dry lines around the eyes.
4. Prevention and control of children's pseudomyopia (eye fatigue).
5. Hydrogen eye massager is an important supplement to dry eye SPA in hospitals. You can also enjoy the wonderful feeling brought by Chinese medicine smoked eyes at home.

For people:

This product is suitable for young students, people who use electronic devices such as mobile phones and computers, and patients with dry eyes and visual fatigue, including students, white-collar workers, overtime workers, elderly people and other family members.


Operation Process Card

①

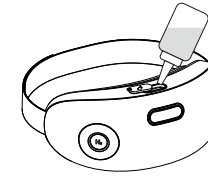


Add mineral water or special atomized eye drops; it is strictly forbidden to use unknown liquids, otherwise the effect will be at your own risk. (It is not recommended to use pure water or distilled water to avoid a small amount of atomization).

②

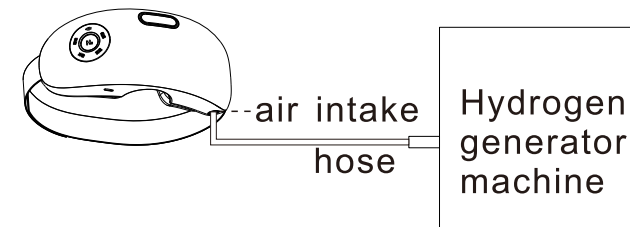

When the gas needs to be implanted, connect the hydrogen generator to the product air inlet with the supplied hose and start it. (The gas inlet end needs to be used with the humidification bottle).

③

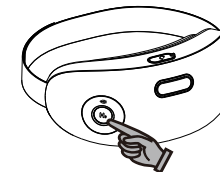

Press and hold the "Switch" button for two seconds to turn on the power. (When "Hot Compress Mode" is selected, please wait for the mist outlet to come out of the mist before wearing it. The warm-up time is about 1-2 minutes).


Maintenance regulations

1. During the warranty period, under normal use of this product, the whole machine will be repaired free of charge for one year (from the date of purchase, the following certificate shall prevail).
2. During the warranty period, if there is one of the following circumstances, the material cost or service fee will be charged as appropriate, please understand;
 - (1) Consumers are damaged due to improper use, maintenance and storage;
 - (2) Damage caused by disassembly by the person who is not responsible for the three-pack repair;
 - (3) Without China's Three Guarantees certificate;
 - (4) The model of the Three Guarantees certificate is inconsistent or altered with the model of the repaired product;
 - (5) The appearance is naturally dirty due to use;
 - (6) Damaged for commercial purposes;
 - (7) Damage caused by force majeure. Not included in the warranty.
3. This product is only suitable for sale and use in mainland China.
4. This card cannot be reissued, please keep it properly.
5. For products that are not within the scope of maintenance and beyond the warranty period, we will still serve you enthusiastically, and only pay for the cost of parts and materials.



Real Steam Experience

Using medical grade ultrasonic atomizer, the size of the mist particles is only 3-5 um, Micron-sized water molecules make the mist more delicate

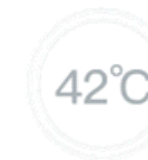

Evenly Steamed at 42°C

Intelligent constant temperature fog heating technology directly heats the fog at a constant temperature, eliminating the bad experience of fluctuating cold and hot, and effectively protecting the eyes.



More Fully Absorbed

The aerodynamic structure is adopted to make the mist output uniformly and at a constant temperature, and directly act on the eyes, so that the eyes are evenly heated, and the mist particles are absorbed more fully.

Product Description



Portable Use

Use Type-c quick charging connector, higher charging efficiency, can be charged and plugged in, easy to carry

Safety Precautions

- ☞ 1. This product is not a substitute for medicine, it can help relieve eye discomfort;
- ☞ Do not take orally;
- ☞ 2. Please be sure to use the atomizing liquid designated by this brand to avoid man-made damage to the instrument due to improper operation;
- ☞ 3. Please follow the operation steps in the manual;
- ☞ 4. Those who are allergic to the components of the atomized liquid are prohibited, and those with allergies should be used with caution;
- ☞ 5. Please follow the doctor's advice when using this product in the treatment of eye diseases;
- ☞ 6. Use the eye spray liquid package within 1 month after opening. If you find turbidity, please stop using it immediately;
- ☞ 7. Do not rinse the main unit of the instrument with water to prevent water from entering the main unit of the product and causing damage to the product;
- ☞ 8. Avoid heavy falls or external impact, and do not disassemble and assemble the host by yourself;
- ☞ 9. Do not use damaged chargers and charging cables to avoid electric leakage, electric shock and fire;

After-sales Service

Warranty Terms

Thank you for choosing this product. In order to better provide high-quality after-sales service, please cooperate as follows:

1. When purchasing, the seller shall issue a sales certificate as the basis for warranty;
2. Please fill in this card in detail and affix the dealer's seal to take effect;
3. Please certify and read the terms of the maintenance regulations.
4. If you need service, please contact your local distributor directly or call the customer service hotline of the headquarters.

Troubleshooting

Unable to heat	Low battery	Recharge
	The atomized liquid in the water tank is used up	Add atomizing liquid
	Incomplete preheating	Wear the product after warming up

Safety Precautions

- ☞ 10. Do not throw this product into fire to prevent explosion;
- ☞ 11. It is forbidden to use this product when its shape changes;
- ☞ 12. User population: Adults, children, children, please ensure that a guardian is present when using it, otherwise incorrect use may cause symptoms to worsen;
- ☞ 13. This product is stored in a well-ventilated room, avoiding high temperature, humidity, and direct sunlight, otherwise it may cause product performance degradation or damage;
- ☞ 14. Prevent infants and people with mental illness from touching when storing or using;
- ☞ 15. If you feel unwell or feel abnormal during use, please stop using it immediately and consult a doctor;
- ☞ 16. Those who have or have had eye diseases, who have had major eye diseases such as surgery, glaucoma, cataracts, retinal detachment, etc., should consult medical staff before using this product;
- ☞ 17. Before use, please remove contact lenses, hairpins, headwear, etc.;
- ☞ 18. If abnormalities or malfunctions are found, please stop using it immediately, disconnect the power supply and contact the after-sales service unit.

Caution


- To avoid the risk of suffocation, please keep children away from the packaging.
- Not applicable to new children and infants.

Safety Precautions

Contraindications

Use Type-c quick charging connector, higher charging efficiency, can be charged and plugged in, easy to carry.

Parts List

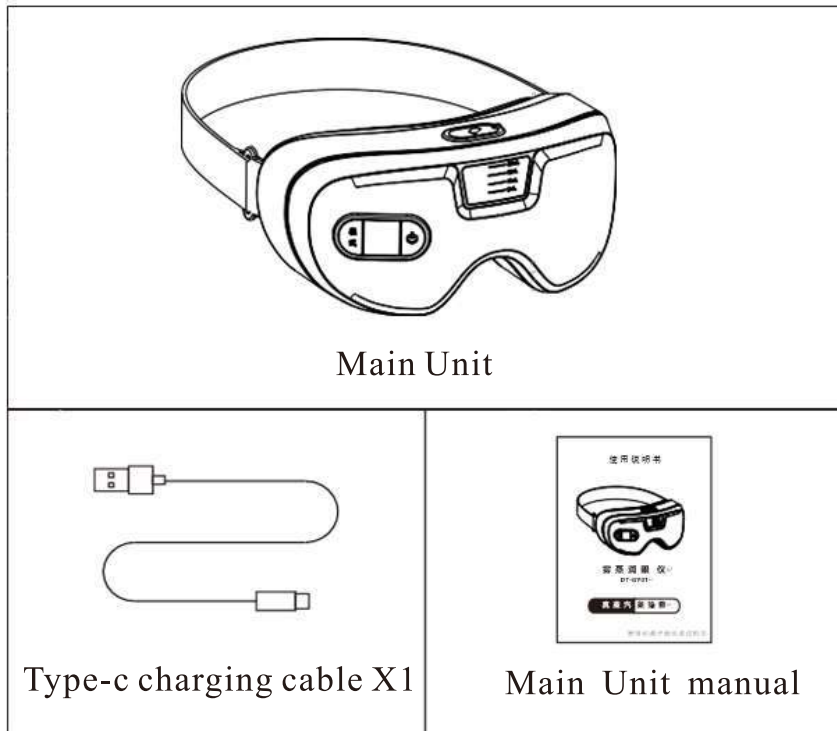

Troubleshooting

Can not boot	Low battery	Recharge
Can't charge	Damaged charging cable	Check or replace the charging cable that meets the product specifications
	Wrong adapter	Check whether the adapter parameters meet the product requirements
Can't get out of the fog	Low battery	Recharge
	Atomizer blockage	Clean the water tank with warm water
	Soft rubber plug vents blocked	Check the vent of the soft rubber plug
	The atomization liquid is too full	Pour out a little atomized liquid
	Blocked bubbles	Shake the product
Eyes have fog on one side and no fog on the other	Condensate blocks the mist outlet	Hold down the mist outlet with your hand, the condensate that blocks the mist outlet will be squeezed out

Maintenance and Storage

Storage

1. When gas needs to be implanted, connect the hydrogen machine (oxygen machine) with the product suction port with a hose
2. In order to avoid damage to the eye instrument, the hydrogen generator (oxygen generator) needs to use the equipment specified by the brand
3. The gas output of the gas Th forming equipment must be less than 1 L/min
4. When the product is in use, please do not block the suction hole.

Tips

- Please place the product out of the reach of children.
- The packaging should be properly stored, and the product can be packed into the packaging when placed.


Troubleshooting

Normal Phenomenon

1. When the product is in use, condensed water will form in the inner cavity of the product, which is normal.
2. When the product is in use, there is a slight vibration, which is normal. Fault self-Check Sheet.

Product Parameter

Name	Hydrogen Eye Massager
Model	MR-03
Product Size	185×65×88mm
Product Weight	360g
Input Voltage	5V
Working Power	<15W
Power Supply	DC3.7V 2000mA
Power Supply	20mL
Ambient Temperature	5°C-40°C
Length Of Work	15min/once
Charging Time	2-3 hours
Executive Standard	GB 4706.1 / GB 4706.15

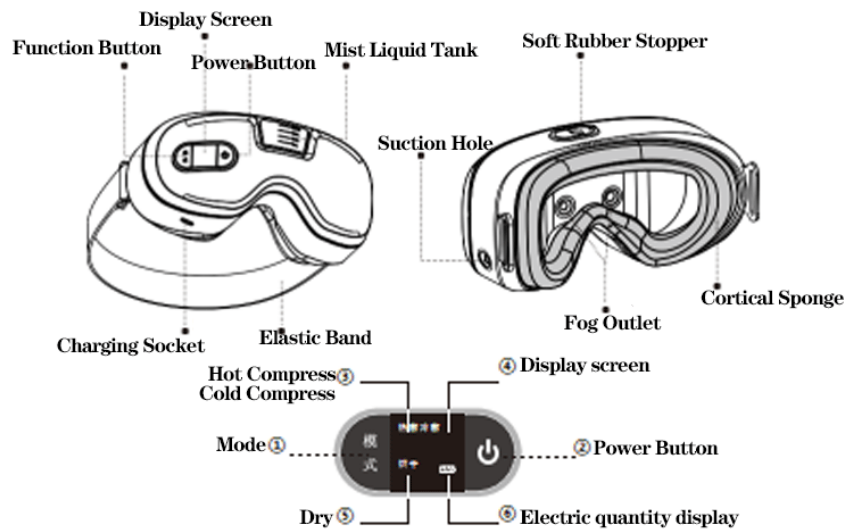

① . Function Selection Buttons

Used to select different function modes, each press changes a function

②. Power Button


Used to turn on/off the power of the product and pause/start when the product is working. Long press this button for 1 second to turn on or turn off the power; when the product is working, press once to pause or start the current working state.

③. Hot Compress Work Light

When this light is on, it indicates that the hot compress function is activated; when it is paused, this light will flash.

Maintenance

Cortical sponge removal method

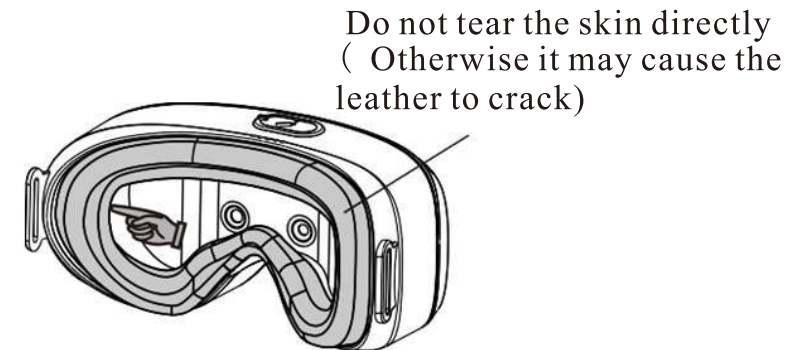

Reach the bottom of the cortex sponge (hard part) with your fingers, and pull out the cortical sponge; when installing, just put the sponge in the corresponding position to stick it.

Use your fingers to reach into the bottom of the leather sponge (hard rubber part), and pull the leather sponge out; when installing, just put the sponge in the corresponding position to stick it.

Maintenance

- ☞ 1. It is forbidden to immerse the product in water or rinse with water.
- ☞ 2. Before cleaning, please turn off the product power and unplug the charging cable.
- ☞ 3. It is forbidden to clean this product with any detergent or cleaning agent.
- ☞ 4. Do not disassemble the product to clean the inside.

Tips

- If the product is dirty, please use a dry soft cloth soaked in clean water and wrung out and wipe the body.
- If necessary, the leather sponge can be removed for replacement or wiping.

④. Cold Compress Work Light

When this light is on, it means that the cold compress function is activated; when it is paused, this light will flash.

⑤. Drying Work Lights


When this light is on, it means the drying function is started; when it is paused, this light will flash. The drying function is used to drain and dry the condensate after the hot compress and cold compress are used.

⑥. Power Indicator

When the power is turned on, the light will be on; when the product is not sufficiently charged, the light will flash rapidly; when charging, the light will flash slowly; when it is fully charged, the light will be on for a long time.

How to use

1. Inject The Atomizing Liquid

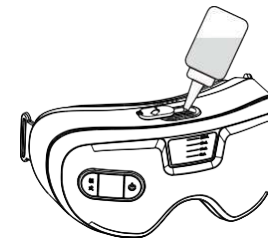

Open the soft rubber stopper and add the atomizing liquid to the mist tank; after adding an appropriate amount of the atomizing liquid, close the soft rubber stopper tightly.

Notes

- In order to avoid bad experience, please use the atomized liquid designated by this brand.
- Do not add tap water or other unknown liquids;
- Each time the atomization liquid is added should not be too full, not more than 20 ml.

2. Turn On The Power



Press and hold the "power button" for 1 second to turn on the power; after turning on the power, the power icon will light up.

3. Function Mode Selection




① Hot Compress Function



1. After the hot compress function is selected, the product has a preheating time of about 30 seconds. During the preheating period, the product does not produce fog and you do not need to wear it immediately; after the preheating is completed, the product emits a warning sound and the fog outlet starts to emit hot fog, and you can wear it now. The product is applied to the eye compresses. (If there is no fog for a long time after the beep sounds, please charge the product before using it).

7. Implant Gas

When this product is used, hydrogen, oxygen and other gases can be implanted according to the characteristics of the atomized liquid or the doctor's suggestion.

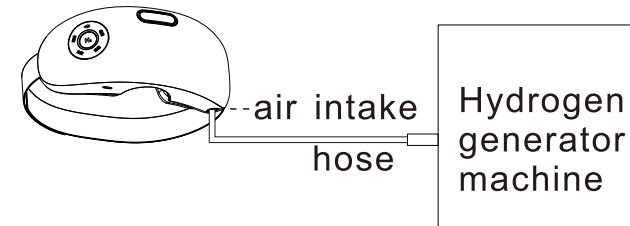

1. When gas needs to be implanted, connect the hydrogen machine (oxygen machine) with the product suction port with a hose
2. In order to avoid damage to the eye instrument, the hydrogen generator (oxygen generator) needs to use the equipment specified by the brand
3. The gas output of the gas forming equipment must be less than 1 L/min
4. When the product is in use, please do not block the suction hole.

6. Charging Instructions

1. Product rated input: 5 V 1 A.
2. Please use a 3C adapter.
3. The adapter should be plugged into a 100-240 V, 50/60 Hz power socket.
4. The charging time is about 2-3 hours, and it can be used 3 times with a single full charge, 15 minutes each time.

Charging



Flashing (charging)



Long on (full charge)

Tips

- Please select a charging adapter that meets the rated voltage of this product.
- Please do not put the product into fire or high temperature place to avoid overheating or explosion of the battery.
- Do not use during charging

① Hot Compress Function



2. The product is set to work for 15 minutes each time, and it will automatically stop working when the time is up.
3. During the hot compress, The condensate will be produced in the inner cavity of the product. After the hot compress is completed, the condensate should be poured out from the mist outlet in time, and the inner cavity of the product should be wiped clean.
4. During the hot compress, if you feel unwell, stop using it immediately.

② Cold Compress Function



1. The product is set to work for 15 minutes each time, and it will automatically stop working when the time is up.
2. During the cold compress process, The condensate will be produced in the inner cavity of the product. After the cold compress is completed, the condensate should be poured out from the mist outlet in time, and the inner cavity of the product should be wiped clean.

How to use


③ Dring Function



- ☞ 1. Drying function: It is used to drain and dry the condensate after hot compress or cold compress.
- ☞ 2. Drying steps: 1. After the hot or cold compress is finished, start the drying function
3. Pour out the condensate with the mist outlet facing down
- ☞ 4. After the condensate is poured out, the product can be placed and allowed to dry automatically, and it will automatically stop after 10 minutes
- ☞ 5. The drying function can effectively keep the inner cavity of the product clean. If you do not use the product for a long time, you should empty the atomized liquid in the mist tank and operate the drying function before placing the product.

4. Wear

Wear Adjust the elastic band to a suitable position and wear the product.

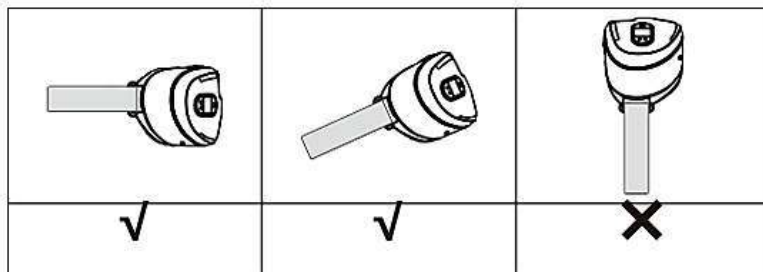

How to use

5. Usage and Dosage

- ①. It is recommended to use 1-2 times a day, put the atomized liquid into the mist liquid tank, turn on the switch and press the function key to use;
- ②. If symptoms such as visual fatigue, dry eyes, and eye pain are sprayed on the eyes, massage the eyes with your hands for 3-5 minutes. The frequency of use can be adjusted appropriately according to the symptoms.
- ③. A single use of 5 g of atomizing liquid can be atomized for 15 minutes, and the amount of liquid filled in the water tank can be used continuously for 4 times. Each bottle of 50 g can be used about 10 times.
- ④. The shelf life of the atomized liquid is 2 years, please use it within 1 month after opening the bottle cap of the atomized liquid;
- ⑤. If the atomized liquid is not used up for more than 3 days after being injected into the water tank, discard the residual eye fluid in the water tank and clean it before use.

Tips

- The product can be paused by pressing the "power button" while it is working. For example, long press the "power button" to shut down.
- When the product is working, please do not turn the product upside down.